

# Status

# DS Nicotine

Rapid immunoassay for the qualitative detection of cotinine in urine.

## Feature


- Simple test procedure
- 500 ng/mL cutoff level
- Detects active nicotine intake
- Helps to permit counseling and intervention immediately



## Test Procedure\*

1. Dispense 3 drops of urine into cassette sample well (S).
2. Read the results after 5 minutes.
3. Negative results: reddish lines appear in both the Control (C) and Test (T) areas.
4. Positive results: Only one colored line appears in the Control area (C).

\*See package insert for details

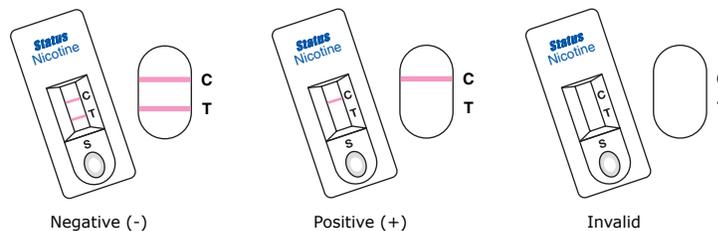

ORDERING INFORMATION		
PRODUCT	# OF TESTS	ITEM #
Status DS Nicotine™	35 tests	21735

For more information or to place an order, contact:



A PBM Group Company

71 Veronica Avenue • Somerset • New Jersey • 08873  
 800.526.2125 • 732.246.3366 • Fax: 732.246.0570 • [www.lifesignmed.com](http://www.lifesignmed.com)